



## Scientific Presentation Integrated BI & CI Technology For The Human Health



LISTER BIOMEDICAL CO., LTD.



# Steam Sterilization

LISTER currently offer three major solutions to ensure that the product is effectively sterilized of steam sterilization: Pre-control, Loads control and Pack control.

## ● Pre-Control

Performed after a warm-up cycle at the beginning of every day, the vacuum-assisted steam sterilizer is in operation, the Pre-control Test verifies the vacuum-assisted sterilizer effectively removes air.

### 1. Type 2 Bowie Dick Test Pack



Non-luminal Units  
**LBD10**

Units: 100pcs/ctn


## ● Loads Control

Monitoring Tool	Frequency
Type 5 CI + PCD	Every batch the sterilizer is run
BI+PCD	Every day the sterilizer is used
BI + Type 5 CI + PCD	For implant load or emergency use load

### 1. Type 5 CI + PCD




Non-luminal Units  
**LPCD12**  
Units: 100packs/ctn

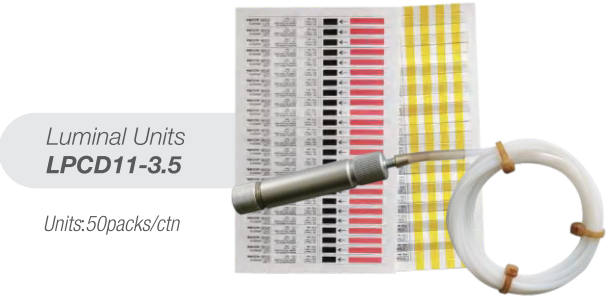

Luminal Units  
**LPCD11-3.5**

Units: 50packs/ctn

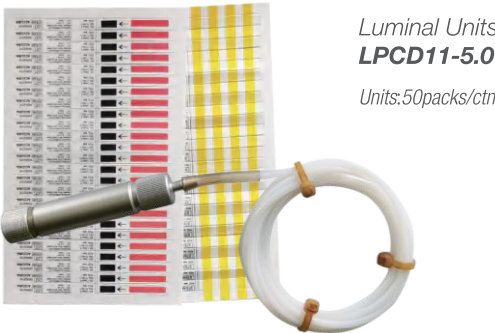

Luminal Units  
**LPCD11-5.0**  
Units: 50packs/ctn

2. BI+PCD



Non-luminal Units  
**LPCD23/24/25**  
Units: 100packs/ctn

3. BI + Type 5 CI + PCD



Non-luminal Units  
**LPCD20/21/22**  
Units: 100packs/ctn

4. BI



**LBS020**  
Units: 50pcs/box  
Chemical Indicator  
Unsterilized  
Sterilized



**LBS060**



**LBS180**

NOTE: If you decide to produce PCD devices that meet the standards yourself, you can also purchase BI and Type 5 CI separately.

Pack Control

The LISTER products listed below are used for Pack-control:

1. Type 1 CI




External  
**LST1250/1950/2550**  
Units: LST1250: 80 Rolls/ctn;  
LST1950: 50 Rolls/ctn;  
LST2550: 40 Rolls/ctn




2. Type 4 CI




Internal  
**LCS40**  
Units: 200pcs/pack

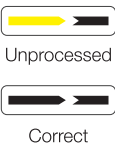

3. Type 5 CI

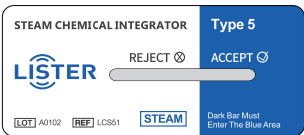

Internal & Batch Control  
**LCS51**

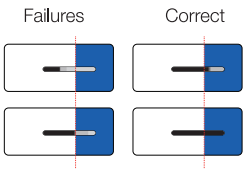

Units: 500pcs/pack



# EO Sterilization

LISTER currently offer two major solutions to ensure that the product is effectively sterilized of EO sterilization: Loads control and Pack control.

## ● Loads Control

Batch control scheme includes biological indicators (BI), chemical indicators (Type 5 CI), and process challenge devices (PCD).  
The LISTER products listed below are used for Batch-control:

### 1. BI



**LBE180**

Units:50pcs/box

Chemical Indicator

Unsterilized

Sterilized

NOTE: A PCD device of certain difficulty is currently under development.A type 5 CI for EO sterilization is currently under development.

## ● Pack Control

The LISTER products listed below are used for Pack-control:

### 2. Type 1 CI



External  
**LET1250/1950/2550**

Units:LET1250:80 Rolls/ctn;  
LET1950:50 Rolls/ctn;  
LET2550:40 Rolls/ctn




Unsterilized



Sterilized

### 3. Type 4 CI

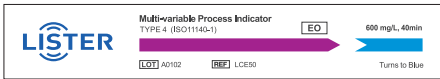

Internal  
**LCE50**

Units:200pcs/pack


Unprocessed



Failures



Correct

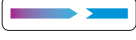



## VH<sub>2</sub>O<sub>2</sub> Sterilization

LISTER currently offer two major solutions to ensure that the product is effectively sterilized of VH<sub>2</sub>O<sub>2</sub> sterilization: Loads control and Pack control.



### ● Loads Control

Batch control for VH<sub>2</sub>O<sub>2</sub> sterilization is the process by which the load is monitored and released based on the result of a Biological Indicator (BI) in a process challenge device (PCD). A PCD contains BI should be used every batch the sterilizer is run, for routine sterilizer efficacy testing.

#### 1. BI



### ● Pack Control

The LISTER products listed below are used for Pack-control:

#### 2. Type 1 CI



External  
**LHT1250/1950/2550**

Units: LHT1250: 80 Rolls/ctn;  
LHT1950: 50 Rolls/ctn;  
LHT2550: 40 Rolls/ctn




Unsterilized



Sterilized

#### 3. Type 4 CI

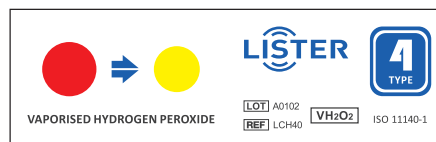

Internal  
**LCH40**

Units: 200 pcs/pack

Unprocessed



Failures



Correct





# FORM Sterilization

LISTER currently offer two major solutions to ensure that the product is effectively sterilized of FORM sterilization: Loads control and Pack control.

## ● Pack Control

The LISTER products listed below are used for Pack-control:

### 2. Type 4 CI

Multi-variable Process Indicator

FORM

LOT

A0102

REF

LCF40

ISO 11140-1

4

TYPE

## ● Loads Control

Batch control for FORM sterilization is the process by which the load is monitored and released based on the result of a Biological Indicator (BI) in a process challenge device (PCD). A PCD contains BI should be used every batch the sterilizer is run, for routine sterilizer efficacy testing.

### 1. BI

Units:50pcs/box

Chemical Indicator

Unsterilized

Sterilized

LBF020

LBF060

Unprocessed


Failures

Correct


Internal  
**LCF40**  
Units:200pcs/pack

# Automatic Biological Reader

The working principle of an automatic biological reader is based on the following chemical reaction:

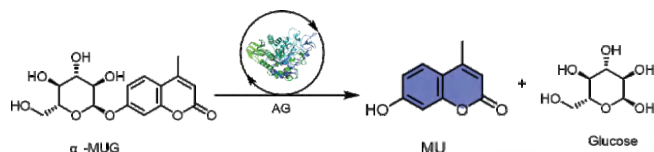

Firstly, after the spore is revived, it will immediately search for an energy source (sugar). In the natural environment, sugar exists in complex polysaccharide forms that microorganisms cannot transport and absorb. Enzyme AG is a powerful assistant for spores, helping them cut complex polysaccharides into monosaccharides.  $\alpha$ -MUG is a disaccharide composed of glucose and fluorescent substance MU. The recovery culture medium in the BI glass tube contains an appropriate amount of  $\alpha$ -MUG.

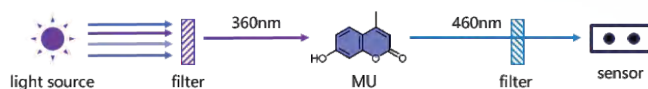

When the culture medium is in full contact with the bacterial carrier, the surviving spores will command AG to cut  $\alpha$ -MUG, and release MU. Individual MU will emit 460nm fluorescence under 360nm light irradiation. Automatic biological reader will determine whether the spores are alive based on the received fluorescence intensity.

The LISTER Automatic Biological Reader listed below:







**BIOPT2**



**BIOPT4**

## Feature Description:

-  Scan the code to obtain information, More Convenient
-  Mobile app remote control, More Freedom
-  Scan the code to obtain information, More Convenient
-  LISTER 2.0 Reading Program, More Accurate



[www.listerbiomed.com](http://www.listerbiomed.com)

